

Escape-Sheet® Standard / Premium

Evacuation directly ... from a bed

What`s an Escape-Sheet®?

People who aren't able to leave their bed independently need to be evacuated using evacuation equipment. One of these aids is the Escape-Sheet®. A mattress with an Escape-Sheet® underneath can be transformed very quickly into an evacuation mattress.

How does an Escape-Sheet® work?

In case of emergency the strong safety belts are taken from underneath the mattress. The belts are strapped over the patient and bedding so the patient is safely secured to the mattress. The Escape-Sheet® can then be easily pulled to the floor and taken to a place of safety.


Advantages Escape-Sheet®:

- Straps and handles highly visible in bright green colour.
- Extra-long towing straps to help operators with improved handling practices.
- Plastic buckle pockets increase the life expectancy of the product (allocated on the side of the sheet).
- Handgrips to assist with sideways manoeuvrability.
- Static-Control™ (anti-static polyester).
- Adjusting corner straps (suitable for each mattress height).
- Premium model features, body-covers and an extra strap at the foot end of the sheet so as to help with easier handling.
- Hydro-tech Nano technology (water-repellent)
- CE & TÜV certificates.

Our advice:

Make sure every bed in your Organization is equipped with an Escape-Sheet®, so you're always prepared for optimal and rapid evacuation.



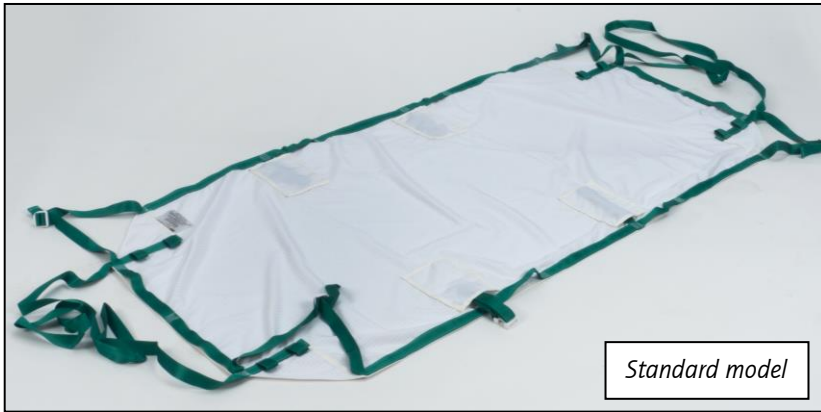

Standard model



Premium model



Extra-long towing straps



Corner sheet adjusters



Buckles in pockets



Option: Straps in 3 different colors

Specifications	Escape-Sheet® Standard	CE	Escape-Sheet® Premium	CE	Escape-Sheet® Custom made	CE
Article number:	21240		21250		21260	
Dimensions (length x width):	± 200 x 90 cm		± 200 x 90 cm		Available on request	
Material:	99% Polyester / 1% Static- Control™		99% Polyester / 1% Static- Control™		99% Polyester / 1% Static- Control™	
Maximum payload:	250 kg		250 kg		250 kg	
Weight:	ca. 1000 gram		ca. 1400 gram		Available on request	
Demonstration and advise provided:	✓		✓		✓	
Ready within seconds:	✓		✓		✓	
Immediate evacuation directly from the bed:	✓		✓		✓	
Industrial wash 75° C:	✓		✓		✓	
Hydro-tech nano technology (water-repellent):	✓		✓		✓	
Adjustable corner straps:	✓		✓		✓	
For all mattress thicknesses:	< 23 cm		< 23 cm		✓	
Body-covers and an extra foot end strap:	-		✓		Available on request	
Buckles in pockets:	✓		-		Available on request	
Heat and flame-resistant version (HFR):	Optional		Optional		Optional	
Towing straps available in 3 different colors (intuitive handling):	Optional		Optional		Optional	
Warranty:	2 years		2 years		2 years	



Interim changes to technical specifications and printing errors reserved.
 Produced by Escape Mobility Production & Innovation: TÜV approved organization..
 Health guaranteed through Öko-tex Standard 100

